

## Product: ZAP70, Active

Catalog #: 02-2062

Amount: 5 µg

### Product Description

Recombinant full-length human ZAP70 was expressed by baculovirus in Sf9 insect cells using a N-terminal GST tag. The gene accession number is [NM\\_001079](#).

### Gene Aliases

SRK; STD; ZAP-70

### Formulation

Recombinant protein stored in 50mM Tris-HCl, pH 7.5, 150mM NaCl, 0.25mM DTT, 0.1mM EGTA, 0.1mM EDTA, 0.1mM PMSF, 25% glycerol.

### Storage and Stability

Store product at -70°C. For optimal storage, aliquot target into smaller quantities after centrifugation and store at recommended temperature. For most favorable performance, avoid repeated handling and multiple freeze/thaw cycles.


### Scientific Background

ZAP70 is a non-receptor protein tyrosine kinase (part of the Syk/Zap70 family) that is involved in signaling by the T-cell antigen receptor (TCR). Ligation of the TCR/CD3 receptor in Jurkat T-cells induces phosphoprotein complexes which contain ZAP70 (1). TCR zeta chains are initially phosphorylated by p56Lck that lead to the recruitment of ZAP70 via its SH2 domain. ZAP70 in turn phosphorylates other proteins in the TCR-phosphoprotein complex. One of the natural substrates for ZAP70 is the zeta-chain dimer of the TCR/CD3 complex (2).

### References


1. Duplay, P. et al: p56lck interacts via its src homology 2 domain with the ZAP-70 kinase. J Exp Med. 1994 Apr 1;179(4):1163-72.
2. Chan, A C. et al: Differential expression of ZAP-70 and Syk protein tyrosine kinases, and the role of this family of protein tyrosine kinases in TCR signaling. J Immunol. 1994 May 15;152(10):4758-66.

### Specific Activity

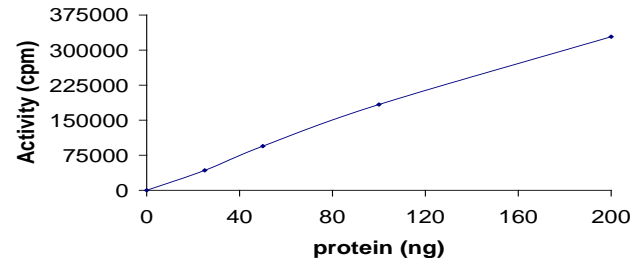

The specific activity of ZAP70 was determined to be **113 nmol /min/mg** as per activity assay protocol.

### Purity

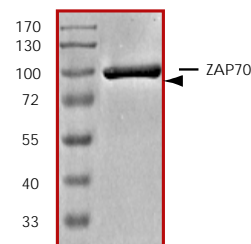

The purity was determined to be **>95%** by densitometry. Approx. MW **96kDa**.

## ZAP70, Active

Full-length recombinant protein expressed in Sf9 cells

Catalog Number	02-2062
Quantity	5µg
Specific Activity	113 nmol/min/mg
Specific Lot Number	B108-1
Purity	>95%
Format	5µg in 50µl
Concentration	0.1µg/µl
Stability	1yr At -70°C from date of shipment
Storage & Shipping	Store product at -70°C. For optimal storage, aliquot target into smaller quantities after centrifugation and store at recommended temperature. For most favorable performance, avoid repeated handling and multiple freeze/thaw cycles. Product shipped on dry ice.

This product is sold for laboratory research use or further manufacturing only and should not be used for human therapeutic or diagnostic applications. The information presented is believed to be accurate; however, said information and products are offered without warranty or guarantee since the ultimate conditions of use and the variability of the materials treated are beyond our control. Nothing disclosed herein is to be construed as a recommendation to use our products in violation of any patents. Under no circumstances shall ARP American Research Products, Inc. be liable for damages, whether consequential, compensatory, incidental or special, strict liability or negligence, breach of warranty or any other theory arising out of the use of the products available from ARP American Research Products, Inc. Nothing contained herein warrants that the use of the products will not infringe on the claims of any patents covering the product itself or the use thereof in combination with other products or in the operation of any process.

## Activity Assay Protocol

### Reaction Components

#### Active Kinase (Catalog #: 02-2062)

Active ZAP70 (0.1µg/µl) diluted with Kinase Dilution Buffer (Catalog #: K21-09) and assayed as outlined in sample activity plot. (Note: these are suggested working dilutions and it is recommended that the researcher perform a serial dilution of Active ZAP70 for optimal results).

#### Kinase Dilution Buffer, pH 7.2 (Catalog #: K21-09)

Kinase Assay Buffer I (Catalog #: K01-09) diluted at a 1:4 ratio (5X dilution) with distilled .

#### Kinase Assay Buffer I, pH 7.2 (Catalog #: K01-09)

Buffer components: 25mM MOPS, 12.5mM β-glycerol-phosphate, 25mM , 5mM EGTA, 2mM EDTA. Add 0.25mM DTT to Kinase Assay Buffer prior to use.

#### [γ]-ATP Assay Cocktail

Prepare 250µM [γ]-ATP Assay Cocktail in a designated radioactive working area by adding the following components: 150µl of 10mM ATP Stock Solution (Catalog #: A50-09), 100µl [32P]-ATP (1mCi/100µl), 5.75ml of Kinase Assay Buffer (Catalog #: K01-09). Store 1ml aliquots at -20°C.

#### 10mM ATP Stock Solution (Catalog #: A50-09)

Prepare ATP stock solution by dissolving 55mg of ATP in 10ml of Kinase Assay Buffer (Catalog #: K01-09). Store 200µl aliquots at -20°C.

#### Substrate

Poly (Glu:Tyr, 4:1) synthetic peptide substrate diluted in distilled to a final concentration of 1mg/ml.

### Assay Protocol

- Step 1. Thaw [32P]-ATP Assay Cocktail in shielded container in a designated radioactive working area.
- Step 2. Thaw the Active ZAP70, Kinase Assay Buffer, Substrate and Enzyme Dilution Buffer on ice.
- Step 3. In a pre-cooled microfuge tube, add the following reaction components bringing the initial reaction volume up to 20µl:
  - Component 1. 10µl of diluted Active ZAP70 (Catalog # 02-2062)
  - Component 2. 10µl of 1mg/ml stock solution of substrate
- Step 4. Set up the blank control as outlined in step 3, excluding the addition of the substrate. Replace the substrate with an equal volume of distilled H2O.
- Step 5. Initiate the reaction by the addition of 5µl [32P]-ATP Assay Cocktail bringing the final volume up to 25µl and incubate the mixture in a water bath at 30°C for 15 minutes.
- Step 6. After the 15 minute incubation period, terminate the reaction by spotting 20µl of the reaction mixture onto individual pre-cut strips of phosphocellulose P81 paper.
- Step 7. Air dry the pre-cut P81 strip and sequentially wash in a 1% phosphoric acid solution (dilute 10ml of phosphoric acid and make a 1L solution with distilled H2O) with constant gentle stirring. It is recommended that the strips be washed a total of 3 intervals for approximately 10 minutes each.
- Step 8. Count the radioactivity on the P81 paper in the presence of scintillation fluid in a scintillation counter.
- Step 9. Determine the corrected cpm by removing the blank control value (see Step 4) for each sample and calculate the kinase specific activity as outlined below.

#### Calculation of [P<sup>32</sup>]-ATP Specific Activity (SA) (cpm/pmol)

Specific activity (SA) = cpm for 5µl [32P]-ATP / pmoles of ATP (in 5µl of a 250µM ATP stock solution, i.e., 1250 pmoles)

#### Kinase Specific Activity (SA) (pmol/min/µg or nmol/min/mg)

Corrected cpm from reaction / [(SA of 32P-ATP in cpm/pmol)\*(Reaction time in min)\*(Enzyme amount in µg or mg)]\*[(Reaction Volume) / (Spot Volume)]

This product is sold for laboratory research use or further manufacturing only and should not be used for human therapeutic or diagnostic applications. The information presented is believed to be accurate; however, said information and products are offered without warranty or guarantee since the ultimate conditions of use and the variability of the materials treated are beyond our control. Nothing disclosed herein is to be construed as a recommendation to use our products in violation of any patents. Under no circumstances shall ARP American Research Products, Inc. be liable for damages, whether consequential, compensatory, incidental or special, strict liability or negligence, breach of warranty or any other theory arising out of the use of the products available from ARP American Research Products, Inc. Nothing contained herein warrants that the use of the products will not infringe on the claims of any patents covering the product itself or the use thereof in combination with other products or in the operation of any process.